

# ASPIRATION CATHETERS

# ASPIRATION CATHETERS

## Export Advance™

### GENERAL CHARACTERISTICS


- Supplied sterile
- Items per box: 1
- Includes two syringes, one filter cup and extension line with stopcock.

### PRODUCT DESCRIPTION

The Export Advance™ aspiration catheter features a preloaded stylet and manual locking syringe to remove embolic material throughout the arterial system.

### FEATURES

- Full Wall Variable Braiding Technology shaft construction
- Preloaded stylet
- Buddy wire compatibility
- Large extraction lumen (0.044" proximal; 0.043" distal)
- Soft, short, forward-facing tip

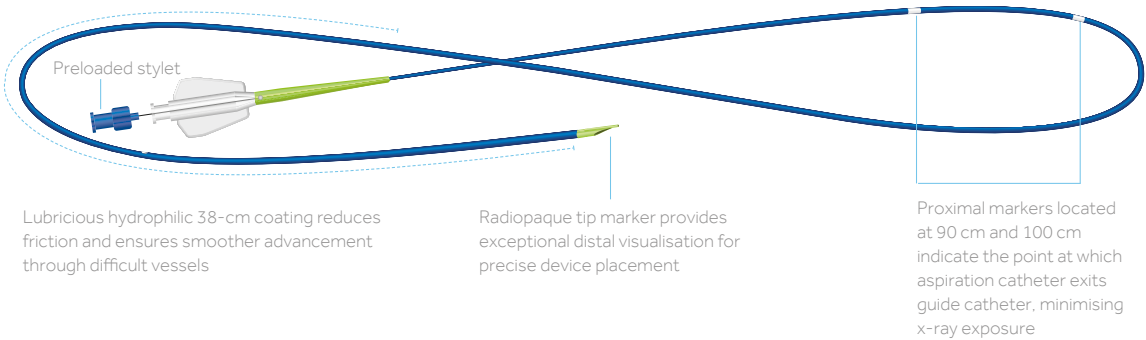

### TECHNICAL INFORMATION

Distal O.D. (in.)	0.067
Lumen I.D. (in.)	0.043 (distal), 0.044 (proximal)
Length (cm)	140
Braiding	Variable, hub to tip
Tip design	Soft, short, beveled tip
Coating	Hydrophilic (38 cm)
Syringe (cc)	2 x 30
Product code	ADVANCECE

# ASPIRATION CATHETERS

## Export™ AP

### GENERAL CHARACTERISTICS


- Supplied sterile
- Items per box: 1
- Includes two syringes, one filter cup and extension line with stopcock.

### PRODUCT DESCRIPTION

The Export™ AP aspiration catheter uses a manual locking syringe to remove embolic material throughout the arterial system.

### FEATURES

- Full Wall Variable Braiding Technology shaft construction
- Extended hydrophilic coating (40 cm)
- Forward-facing atraumatic tip

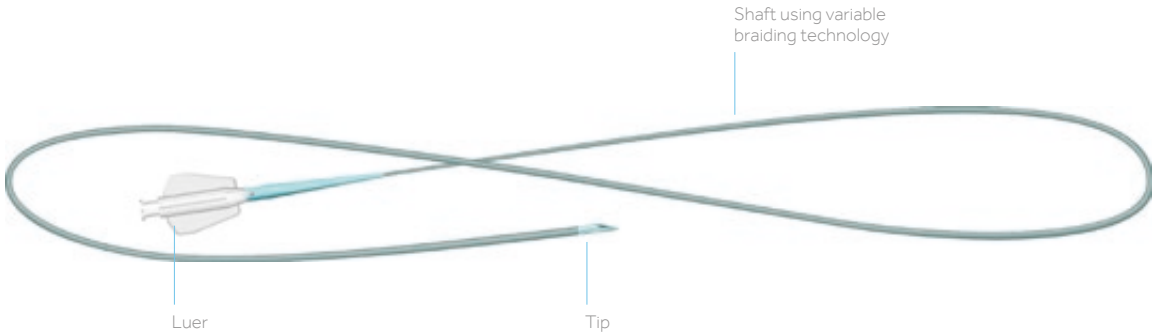

### TECHNICAL INFORMATION

Lumen I.D. (in.)	0.043
Lumen O.D. (in.)	0.068
Guide catheter compatibility (in.)	6F (min ID 0.070)
Guide wire compatibility (in.)	0.014
Delivery System	Rapid Exchange
Length (cm)	140
Radio-opaque marker	1.5 mm from distal tip
Product code	EXPORTAPCE

# ASPIRATION CATHETERS

## Export™ 7 F

### GENERAL CHARACTERISTICS

- Supplied sterile.
- Includes two syringes and two filters.

### PRODUCT DESCRIPTION

The Export™ Aspiration Catheter uses a manual locking syringe to remove embolic material throughout the arterial system. Export™ 7 French delivers double the aspiration power.

### FEATURES

- Distal lumen (in.) 0.050
- Lumen O.D. (in) 0.078
- Guide catheter compatibility (in.) 7F (min ID 0.080)



### TECHNICAL INFORMATION

PRODUCT CODE	GUIDE COMPATIBILITY (IN.)	WIRE COMPATIBILITY (IN.)	LENGTH (CM)
GEZ62007B	7F min. guide I.D. 0.080	0.014	145